

The Classic-Check® System

Primary applications: Beds and bedside chairs (AC powered)

Through thirty years of designing and manufacturing fall and restraint reduction monitoring systems, Bed-Check® has learned what caregivers want most: simplicity, versatility, and reliability. Classic-Check® has classic features for monitoring fall risk patients and residents:

Day/Night Mode – Volume is pre-programmable to different volumes for Day and Night. Press one button to change modes.

Night Light – Soft, blue light within the splash guard (when in Night Mode).

Easily distinguished from other alarms:

Ten programmable alarm tones – Six are tunes, which may be more suitable for easily agitated patients.

Nurse Call Capability – Connect to virtually any known nurse call system (meets UL® 1069). One version comes with a port for the nurse call button. The other, without that port, requires a “Y” adaptor to share a call system wall port with the call button.

For Hospital and LTC – Utilizes an IT (Intermediate Term) or LT (Long Term) Sensormat®.

“Just right” volume – Seven settings, 63-85db, to suit the distance from patient/resident to caregiver, and separately settable for Day and Night Modes.


Servicing the patient/resident – Put the System on hold for 25 seconds while servicing or removing the patient. Resumption of monitoring is automatic.

Set delay to suit your patient/resident – Settable for 1-3 seconds. Longer for restless individuals (reducing “false” alarms), and shorter for fast exiting patients.

Mounting choices – Velcro® mounting strap, or optional wire holster for headboard or wall mounting.

Static protection – Designed/tested to withstand high environmental static voltage commonly encountered in healthcare surroundings (>22kV).

Integral cords – Integral 10’ power and nurse call system cords. (Also available without nurse call cord.)

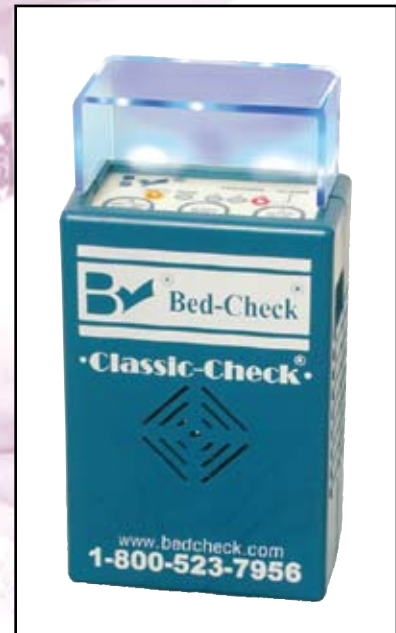

Classic-Check® Control Unit

